

Iron Assay Kit (Colorimetric)

LS-K208-250 (250 Tests) • See Storage Conditions Below



Introduction

Iron level in blood is a reliable diagnostic indicator of various disease states. Increased levels of iron concentration in blood are associated with blood loss, increased destruction of red blood cells (e.g. hemorrhage) or decreased blood cell survival, acute hepatitis, certain sideroachrestic anemias, ingestion of iron-rich diets, defects in iron storage (e.g. pernicious anemia). Decreased levels of blood iron may result from insufficient iron ingestion from diets, chronic blood loss pathologies, or increased demand on iron storage as during normal pregnancy.

Simple, direct and automation-ready procedures for measuring iron concentrations find wide applications in research, drug discovery and environmental monitoring. This iron assay kit is designed to measure total iron directly in serum without any pretreatment. The improved Method utilizes a chromogen that forms a blue colored complex specifically with Fe^{2+} . Fe^{3+} in the sample is reduced to Fe^{2+} , thus allowing the assay for total iron concentration. The intensity of the color, measured at 590nm, is directly proportional to the iron concentration in the sample.

Key Features

- Sensitive and accurate. Linear detection range 27 $\mu\text{g}/\text{dL}$ (4.8 μM) to 1,000 $\mu\text{g}/\text{dL}$ (179 μM) iron in 96-well plate assay.
- Simple and high-throughput. The procedure involves addition of a single working reagent and incubation for 40 min. Can be readily automated as a high-throughput assay for thousands of samples per day.
- Improved reagent stability and versatility. The optimized formulation has greatly enhanced reagent and signal stability. Cuvette or 96-well plate assay.
- Low interference in biological samples. No pretreatments are needed. Assays can be directly performed on serum samples.

Applications

- Direct Assays: iron in biological samples (e.g. serum).
- Drug Discovery/Pharmacology: effects of drugs on iron metabolism.
- Environmental Monitoring: iron in soil extracts, mineralized samples.

Components

Component	K208-250
	250 Tests
Reagent A	50 mL
Reagent B	4 mL
Reagent C	4 mL
Iron Standard (10 mg/dL Fe^{2+})	1 mL

FOR RESEARCH USE ONLY! Not for use in humans.

LifeSpan BioSciences, Inc. • 2401 Fourth Avenue, Suite 900, Seattle, WA 98121
www.LSBio.com • (206) 464-1554 • TechnicalSupport@LSBio.com

Iron Assay Kit (Colorimetric)

LS-K208-250 (250 Tests) • See Storage Conditions Below



Materials Not Supplied

Pipetting devices and accessories.

Procedure using 96-well plate: Clear bottom 96-well plates (e.g. Corning Costar) and plate reader.

Procedure using cuvette: Cuvettes and spectrophotometer for measuring OD at 510-630nm.

Storage

The kit is shipped at room temperature. Store Reagent A at room temperature and other reagents at 4 °C. Shelf life: 12 months after receipt.

Assay Procedure

Note: (1). Iron chelators (e.g. EDTA) interfere with this assay and should be avoided in sample preparation. (2). Serum or plasma samples should be clear and free of precipitates or turbidity. If not, centrifuge or filter to clarify samples prior to assay. (3). This kit can be applied to measure Fe²⁺ (vs. total iron) content. Prepare Working Reagent by mixing 20 vol of Reagent A, 1 vol of water and 1 vol Reagent C (no reductant in the Working Reagent). The procedure is the same as described for the total iron assay.

Procedure using 96-well plate

- Standards. Prepare 400 µL 1000 µg/dL Premix by mixing 40 µL 10 mg/dL standard and 360 µL distilled water. Dilute standards as follows.

No	Premix + H ₂ O	Vol (µL)	Fe (µg/dL)
1	100µL + 0µL	100	1000
2	80µL + 20µL	100	800
3	60µL + 40µL	100	600
4	40µL + 60µL	100	400
5	30µL + 70µL	100	300
6	20µL + 80µL	100	200
7	10µL + 90µL	100	100
8	0µL + 100µL	100	0

Transfer 50 µL diluted standards and 50 µL sample into a clear flat bottom 96-well plate. For serum/plasma samples, it is recommended to run a sample blank (i.e. a 50 µL sample in a separate well).

- Prepare enough Working Reagent by mixing 20 volumes of Reagent A, 1 volume Reagent B and 1 volume Reagent C. Fresh reconstitution is recommended. Equilibrate to room temperature before assay.

Add 200 µL Working Reagent to Standards and Samples wells. (For serum/plasma samples which require a Sample Blank Control, add 200 µL Reagent A to the Sample Blank wells). Tap plate to mix.

- Incubate 40 min at room temperature and read optical density at 510-630nm (peak absorbance at 590nm).

Procedure using cuvette

- Prepare standards as in 96-well assay. Set up centrifuge tubes labeled Standards and Samples. Transfer 250 µL standards and samples to tubes.
- Add 1000 µL Working Reagent to all tubes. Mix by vortexing. Incubate 40 min at room temperature.
- Transfer to cuvettes and read OD at 590nm (510nm-630nm).

FOR RESEARCH USE ONLY! Not for use in humans.

LifeSpan BioSciences, Inc. • 2401 Fourth Avenue, Suite 900, Seattle, WA 98121

www.LSBio.com • (206) 464-1554 • TechnicalSupport@LSBio.com

Iron Assay Kit (Colorimetric)

LS-K208-250 (250 Tests) • See Storage Conditions Below

Calculations

Subtract OD of “0 µg/dL Fe” from all other standard OD values and plot the OD against standard iron concentrations. Determine the slope using linear regression fitting. Iron concentration of the sample is calculated as


$$[\text{Iron}] = \frac{\text{OD}_{\text{SAMPLE}} - \text{OD}_{\text{BLANK}}}{\text{Slope}} \quad (\mu\text{g/dL})$$

Where OD_{BLANK} is OD values of the water blank (Standard #8), or Sample Blank, if a sample blank is used (e.g. serum or plasma). Typical serum iron values: 70-180 µg/dL.

Conversions: 1 mg/dL Fe equals 179 µM, 0.001% or 10 ppm.

Sample Data

Mouse serum, fetal bovine serum (Invitrogen), and goat serum (Invitrogen) were assayed using the 96-well plate assay protocol. The iron concentrations were 173 ± 2 (n = 4), 149 ± 1 (n = 4), 88 ± 2 µg/dL (n = 4), respectively. Coefficient of variance < 2%.

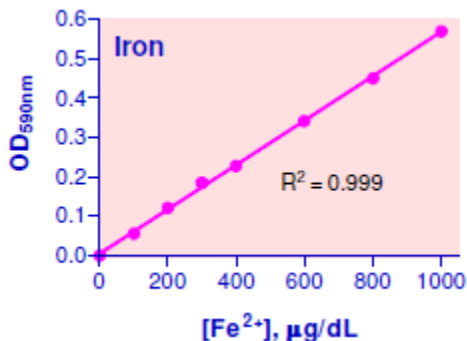

Standard Curve in 96-well plate assay

Version: V.08.09.2018

FOR RESEARCH USE ONLY! Not for use in humans.

LifeSpan BioSciences, Inc. • 2401 Fourth Avenue, Suite 900, Seattle, WA 98121

www.LSBio.com • (206) 464-1554 • TechnicalSupport@LSBio.com